



better analysis counts

Анализатор серы Sindie 262 G3 M

Техническое описание

Архангельск (8182)63-90-72
Астана +7(7172)727-132
Белгород (4722)40-23-64
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89
Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Казань (843)206-01-48

Калининград (4012)72-03-81
Калуга (4842)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Липецк (4742)52-20-81
Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41

Нижний Новгород (831)429-08-12
Новокузнецк (3843)20-46-81
Новосибирск (383)227-86-73
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Пермь (342)205-81-47
Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78

Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Тверь (4822)63-31-35
Томск (3822)98-41-53
Тула (4872)74-02-29
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Уфа (347)229-48-12
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Ярославль (4852)69-52-93

16 32.06

Sindie 2622 G3 M-SERIES

Sulfur Analyzer



Sulfur Analysis with Compliance Flexibility

The Sindie 2622 bench-top analyzer complies with ASTM D2622, D7039 and ISO 20884 methods, enabling complete flexibility in sulfur analysis. No compromises in detection, performance and reliability - the SINDIE 2622 analyzer is the ideal sulfur analytical solution from ultra low sulfur diesel and gasoline to heavy fuel oil and crudes. One analyzer – three compliance solutions.

Application Areas:

- Total sulfur analysis from ultra low sulfur fuels to crudes.
- For use in refinery labs, pipeline terminals, additive plants and inspection laboratories.
- Complies with ASTM D2622, D7039 and ISO 20884.

Features and Benefits:



- LOD: 0.15 ppm at 300 s.
- Dynamic Range: 0.15 ppm to 10%
- Fits on any bench:
37 cm (w) x 50 cm (d) x 34 cm (h).
- Plug-it-in and measure: power is the only utility.
- Touch Screen user interface.
- User programmable measurement time: 30-900 s.
- Two calibrations cover a variety of products over full dynamic range:
 - 0.15 – 3000 ppm wt
 - 0.3 – 10% wt
- No conversion gasses, heating elements, quartz tubes or columns.
- 75 W air-cooled excitation tube.
- Robust polyamide window for easy cleaning.

Options:

- LIMS data output software capability.

MWD XRF

Monochromatic Wavelength Dispersive X-Ray Fluorescence (MWD XRF) utilizes state-of-the-art focusing and monochromating optics to increase excitation intensity and dramatically improve signal-to-background over high power traditional WD XRF instruments. This enables significantly improved detection limits and precision and a reduced sensitivity to matrix effects. A monochromatic and focused primary beam excites the sample and secondary characteristic fluorescence x-rays are emitted from the sample. A second monochromating optic selects the sulfur characteristic x-rays and directs these x-rays to the detector. MWD XRF is a direct measurement technique and does not require consumable gasses or sample conversion.

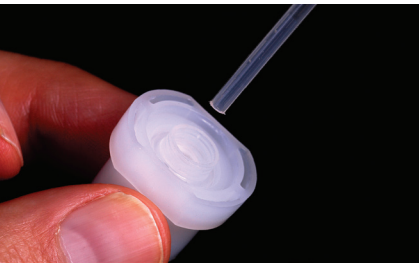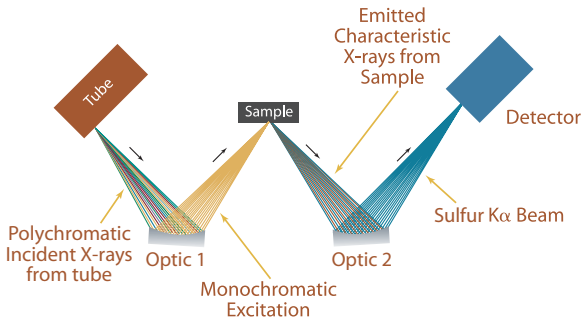

ACCU-CELL Sample Cups


- No assembly of separate film & cup components
- Pre-vented sample cups
- Eliminates sample & cup contamination
- One discharge of 1 ml pipette will fill the cup

Precision


Typical repeatability (r) and reproducibility (R) values in diesel fuel, at 95% confidence. 300 s measurement time.

Sulfur Concentration (ppm)	r	R
2	0.3	0.7
5	0.5	0.8
8	0.6	1.0
15	0.8	1.4
100	2	4
500	5	10

SINDIE-XR Calibration Curve

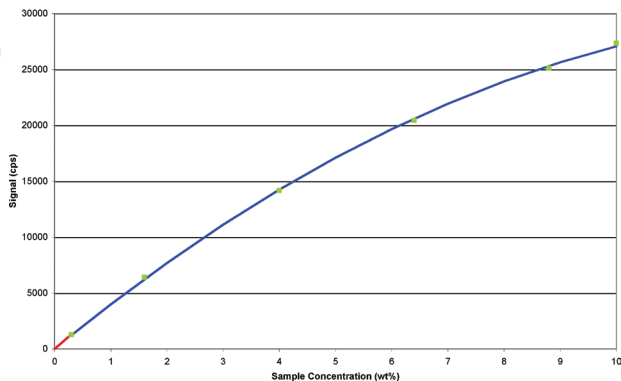

Low Range Calibration

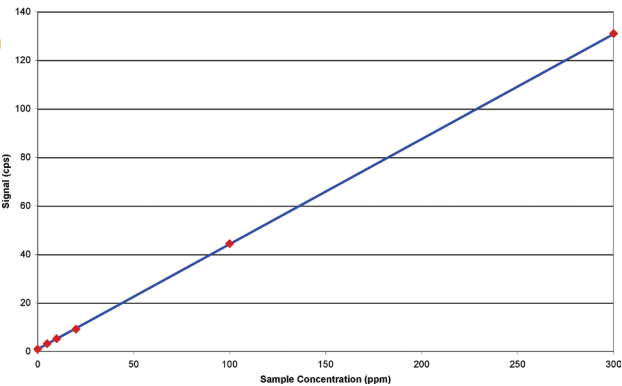

Product Specifications

Test Method	ASTM D7039, D2622 and ISO 20884
Dimensions	37 cm (w) x 50 cm (d) x 34 cm (h)
Power	100-120 VAC, 47-63 HZ at 6.0 Amps/200-240 VAC, 47-63 HZ at 6.0 Amps
Sample Cup Volume	1 ml
I/O Ports	Ethernet 10/100 base T, RS232
Optional Computer Interface	Pentium, 100 MHz, 32 MB RAM/Windows 98 or newer operating system
Ambient Temperature Requirements	5-40° C (40-104° F)
Dynamic Range	Standard: 0.15 ppm – 10% ppm
Measurement	User selectable: 30-900 s
Calibration	8 calibration curves. Automatic and Manual Calibration functionality



better analysis counts

Архангельск (8182)63-90-72
Астана +7(7172)727-132
Белгород (4722)40-23-64
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89
Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Казань (843)206-01-48

Калининград (4012)72-03-81
Калуга (4842)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Липецк (4742)52-20-81
Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41

Нижний Новгород (831)429-08-12
Новокузнецк (3843)20-46-81
Новосибирск (383)227-86-73
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Пермь (342)205-81-47
Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78

Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Тверь (4822)63-31-35
Томск (3822)98-41-53
Тула (4872)74-02-29
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Уфа (347)229-48-12
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Ярославль (4852)69-52-93